



# TRIOS™

Advanced IPL treatments for enhanced customer satisfaction.

Integrating Trios™ into your clinic is the smart decision – for you and your clients. The ability to perform the most widely requested aesthetic treatments – hair removal, skin rejuvenation and acne clearance – gives you a decided edge against your competitors. Trios™ is a leading IPL device that delivers safe and effective outcomes, allowing you to benefit from Viora's extensive expertise in the industry and your patients to enjoy the long-lasting results. By offering a wider spectrum of solutions to your clients, you can increase customer satisfaction and boost revenues.

MK-006 CB1 UK B

### About Viora

Viora, a leading provider of medical aesthetic solutions, is passionate about helping customers deliver the most effective results to their patients. Viora improves upon treatment experiences by enhancing existing technology to give customers the best non-invasive solutions available. With a presence in over 65 countries, Viora's non-invasive products address a wide range of symptoms, treating today's most commonly requested applications for anti-aging, body & facial contouring, cellulite reduction, acne clearance, stretch mark reduction, hair removal, skin lightening and more.


[www.vioramed.com](http://www.vioramed.com) | [info@vioramed.com](mailto:info@vioramed.com) |  [vioramed](#)

**Manufacturer's Address**  
450 Bath Road, Longford  
Heathrow UB7 0EB, UK  
Tel: +44 (0) 208 757 8750  
Fax: +44 (0) 208 757 8767

**Contact Information:**  
**Viora Inc.**  
392 Summit Ave.  
Jersey City, NJ 07306 USA  
Tel: +1-201-332-4100  
Fax: +1-201-332-4200



FDA cleared for marketing  
This device complies with the European Medical Device Directive 93/42/EEC

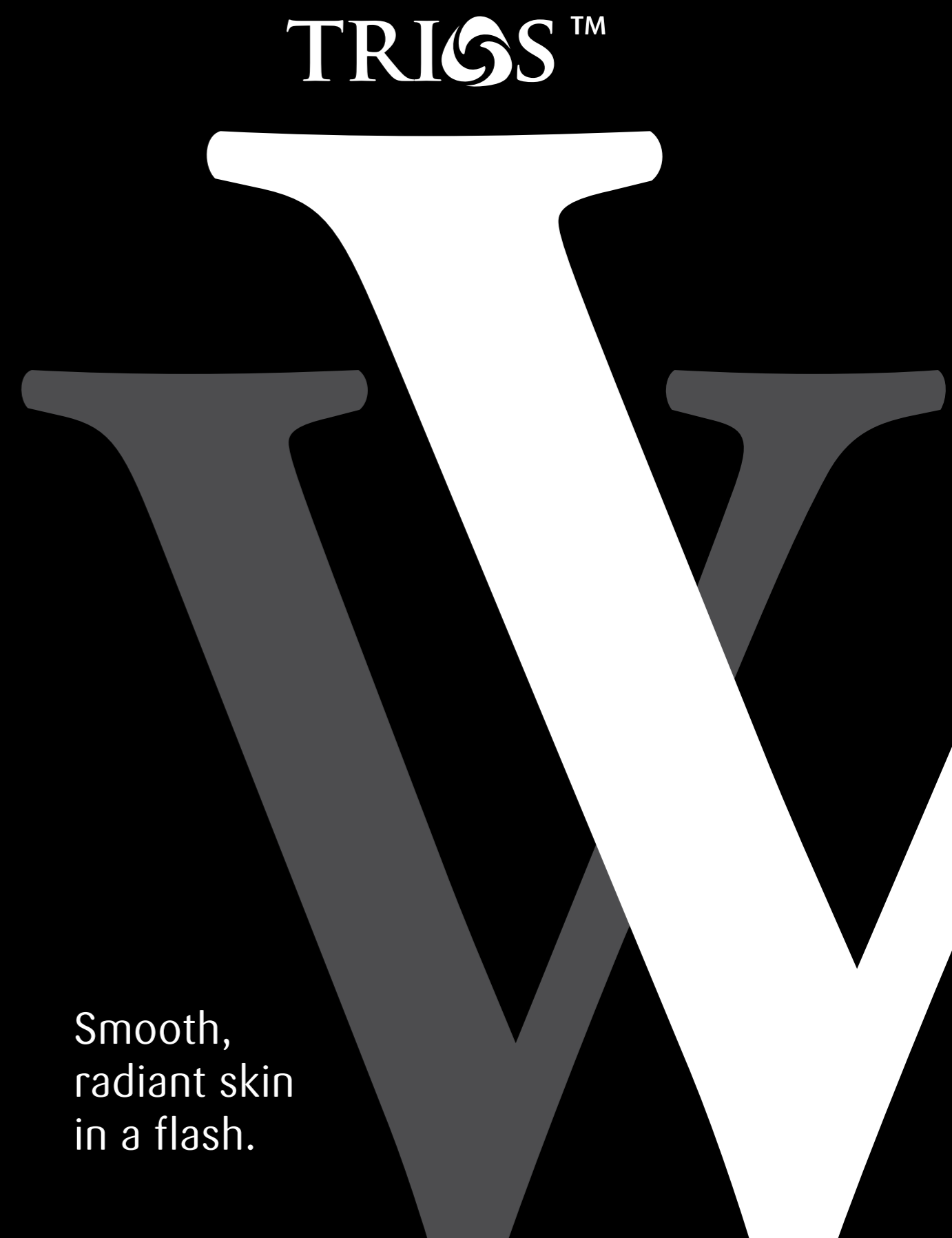

# TRIOS™

Smooth,  
radiant skin  
in a flash.

## The Ability to Offer More

IPL devices are known for being the workhorses of the aesthetic industry. Providing solutions for the most requested applications around the world, such as hair removal, acne clearance and skin rejuvenation, IPL is practically a necessity in every clinic.

Trios™ by Viora, developed by an R&D team with over two decades of engineering experience in IPL based systems, is a modular system with versatile applications that enable you to expand your treatment portfolio and satisfy your customers' demands.

Trios offers a full suite of aesthetic solutions:

- Long-term hair removal for normal hair growth & Hirsutism
- Rejuvenates skin, reducing vascular and pigmented lesions
- Smooths fine lines and wrinkles
- Balances skin pigmentation caused by sun damage & age spots
- Acne clearance



*"The majority of clientele in my business used to be women from 35-50 years of age but since I purchased the Trios, I'm now able to offer so much more to young women, men of all ages and teenagers as well. Trios™ has really upgraded my business and I'm able to give my clients exactly what they need at a great price."*

*- C.W., aesthetician*

## Results, Right on Target

Based on a process called selective photothermolysis, Trios™ emits pulsed light energy in one single flash. The light emitted is absorbed by targeted chromophores (light sensitive cells) and encourages a specific biological process in order to achieve the desired clinical result.



## A Host of Advantages

Trios™ was designed to ensure you are equipped with a system that delivers uncompromising performance, safety and efficiency on all levels.

Trios™ offers advantages that are second to none on the market:

- Actual 22 J/cm<sup>2</sup> of power at skin level for effective treatments and maximum results.
- Highly accurate program settings enable you to tailor the best treatment for your patients.
- Highly reliable device; manufactured to the highest standards and requires minimal setup & maintenance.
- Safety in mind: advanced computerized management control to ensure system performance and risk-free treatments.
- Visible results after the first treatment.
- Fast, comfortable treatments.
- Digital, user-friendly interface.
- Unique, ergonomic hand piece design.
- Portable device allows users to move between rooms and salon locations.

## The Profitable Formula

### Large spot size area (7.5cm<sup>2</sup>)

- Reduced treatment time
- Increased accuracy and efficacy
- Less pulses required per treatment

### Fast repetition rate

- Faster treatment times
- Perform more treatment procedures per day

### Long lamp life (100,000 pulses)

- Cost efficient

### Multiple treatment applications

- Attract new customers
- Increase treatment offering to existing customers

Viora's approach to your success doesn't end with supplying the highest quality systems on the market. We are committed to our customers and continue to pursue clinical, technical and business pathways that will lead to your success.



## Viora Support

Viora offers the best marketing, clinical and technical support available. These ongoing tools and techniques allow you to be secure in the knowledge that a stable, committed company is partnering with you, for your success.

## Technical Specifications

Output	Up to 22 J/cm <sup>2</sup>
Wavelength	HR: 530-1200 nm SR: 570-1200 nm AC: 400-1200 nm
Pulse Duration	25 msec
Pulse Repetition Rate	0.33Hz (every 2-3 seconds)
Spot Size	15 X 50 mm (7.5cm <sup>2</sup> )
Dimensions	45cm (L) x 38cm (W) x 27cm (H)
Weight	17 kg / 37 lbs
Electrical	220V/50Hz or 110V/60Hz - 800 V/A

